Home
Dental Center
RST5010Smart

Please read the product manual carefully
before use.

USER MANUAL

Important safeguard.....	2
Configuration.....	4
Getting started.....	5
How to use Oral irrigator.....	5
How to use Sonic Toothbrush.....	7
How to use Dental Mirror.....	8
Maintenance.....	9

Please read the safety instructions before using this product. The safety instructions provided here help you use this product safely and correctly, and also to prevent you from injuring yourself or others. When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

READ ALL INSTRUCTIONS BEFORE USING

Danger–Water

- Do not place or store the product where it can fall or be pulled into bathtub, sink or toilet.
- Do not place or drop the charger into water or other liquid.
- Do not reach for the charger base that has fallen into water or other liquid. Unplug immediately.
- Do not use while bathing.
- Unplug the charger base before cleaning and ensure that it is dry before plugging it into an electrical outlet.

Danger–Voltage

- Do not plug the charger base into an outlet with a voltage other than specified on the bottom of the charger(100–240V).
- Voltage converters do not guarantee compatibility.
- Do not force the plug into an outlet.
- If the plug does not fit easily into the outlet, discontinue use.
- This appliance has a polarized plug (one blade is wider than the other). As a safety feature, this plug will fit in a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician. Do not attempt to defeat this safety feature.

Medical warning


- If you have had oral surgery or gum surgery within the past two months, consult your dentist before using the product.
- Consult your dentist if excessive bleeding occurs while using this product or minor bleeding continues to occur after one week of use.
- Consult your physician or implant device manufacturer prior to using this product if you have questions concerning the combined use of this product with a pacemaker or other implanted medical device(s).

WARNINGS

To reduce the risk of burns, electrocution, fire, or injury to persons:

- This product is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities or lack of experience and knowledge, unless they have been given adequate supervision or instruction concerning use of the device by an adult responsible for their safety.


- Children should be supervised to ensure they do not play with the product.
- Use this product only for its intended use as described in this manual. This product is designed to clean your teeth, gums and tongue only, Do not use any attachments or brush heads other than those recommended by the manufacturer.
- If your toothpaste includes peroxide, baking soda or other bicarbonate (common in whitening toothpastes), please thoroughly clean your brush head and handle to avoid cracking the plastic.
- Certain liquid mouthwash may damage the device by causing the housing or tank to crack, when using mouthwash, rinse the tank and main unit with water immediately after use. Be especially careful when washing the hose connections at the main unit and the tank.
- Never fill the tank with hot water! The gum lines can be very sensitive and can be easily damaged by hot water.
- Use the unit within an ambient temperature range of 32F–104F
- Never operate this product if it has a damaged cord, plug or charger. Do not use this product if it is not working properly, if it has been dropped or damaged or dropped into water. Immediately discontinue use if the product appears damaged in any way.
- Keep the product away from heated surfaces.
- If you experience discomfort or pain in your mouth, discontinue use of this product and contact your dentist.
- Do not clean the product in the dishwasher.
- Consult your physician prior to using the toothbrush if you have any medical concerns.
- The toothbrush has been tested and is certified as per ETL regulations
- Never drop or insert any object into any opening.
- For household use only. Do not use outdoors.

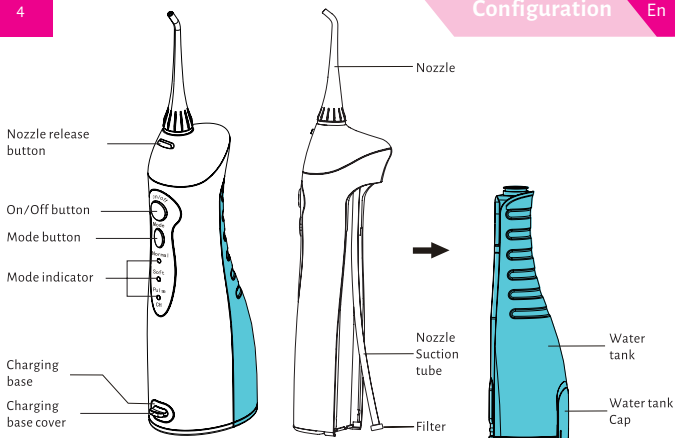
ACCESSORY

BRUSH HANDLE

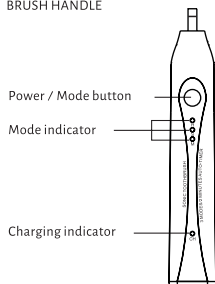CHARGER

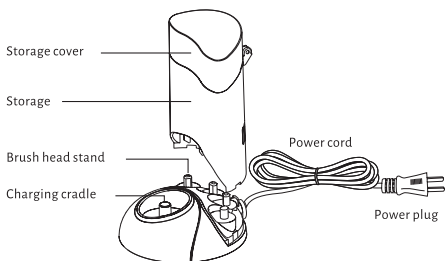
Charging the device(Fig.1)

1. Take the charger base cover down and insert the appliance plug securely into the appliance socket (make sure no liquid remaining inside the socket before insert the appliance plug).
2. Insert the charging adaptor securely into the power outlet.
3. The “CH” indicator blinks when the device is charging and it will turn off after it's fully charged.
4. Recover the cover on to charger base again after charging finished.

Note:


A full recharging will take about 8 hours.

Fig.1

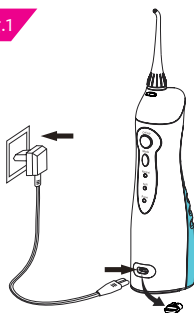How to use Oral irrigator

Fig.2

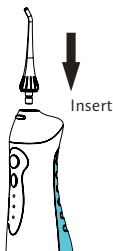1. Attach nozzle onto the unit. (Fig.2)
2. Open the water tank cap to fill the tank with water which temperature not exceed 40°C then close the cap.(Fig.3)
3. Or, Slide the water tank down off the unit, fill the tank with water which temperature not exceed 40°C, then attach the tank onto the unit. (Fig.4)
4. Press “Mode” key to select the mode you desired. During irrigating, you can change the working mode. (Fig.5)

Fig.3


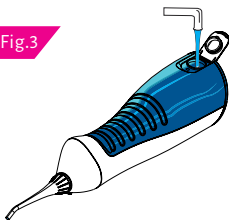
Fig.4

Slide down

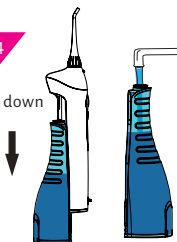Normal: for normal cleaning.

Soft: for gentle cleaning of sensitive gums.

Pulse: for cleaning and massaging the gums.

5. Hold the unit in an upright position with nozzle tip point toward your teeth.

6. Keep your mouth opened slightly so that the water can flow out.

7. Press "On/Off" key to start irrigating. (Fig. 6)

1). Directing the water jet flow at a 90-degree angle to the teeth and gums

2). Directing the water jet slowly along the teeth

3). Directing the water jet along the gum line with the nozzle tip at a 90-degree angle to the gum tissue

4). Cleaning between the teeth and around any braces, crowns or bridgework thoroughly.

8. Press "On/Off" key and then press the nozzle release button to replace other nozzle when you need. (Fig.7)

9. Once irrigating is complete, press "On/Off" key again to turn off the device.

After using the appliance (Fig. 8)

1. Open the water tank cap and dispose of the water in the tank.

2. Press "On/Off" key to turn on the unit, and drain the water inside the appliance.

3. Press "On/Off" key again to turn off the unit, and close the water tank cap.

4. Wipe the device with a dry cloth.

Note: To prevent bacteria propagating, make sure there is no remaining water inside the unit.

If you're not going to use this appliance for a long time, be sure wipe and dry the appliance ASAP before storing it.





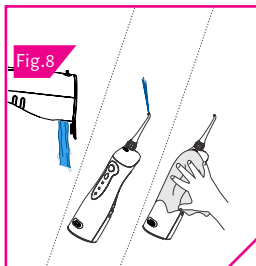
Fig.9

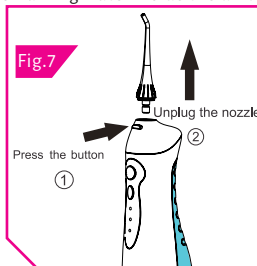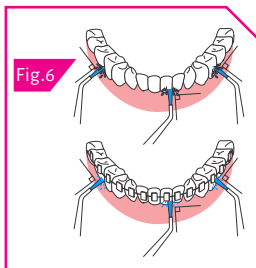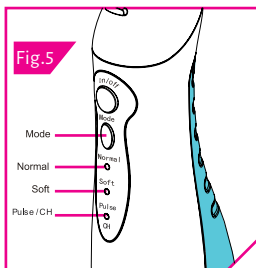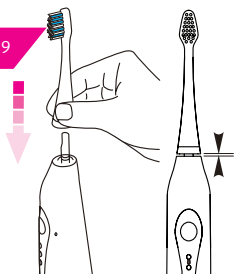1. Slide the brush head onto the handle. Push brush head down until it locks into place on handle.

Notes: There is small gap between the brush head and the handle.(Fig.9)

2. Wet the bristles and apply a small mount of toothpaste.

3. Place the bristles against the teeth at a slight angle towards the gum line.(Fig. 10)

4. Press the Power / Mode button to turn on the toothbrush, It starts working at Normal mode.

Normal(H): for normal cleaning.

Soft(S): for gentle cleaning of sensitive gums.

Massage (M): for cleaning and massaging the gums.

Notes: During brushing, if press the Power / Mode button, toothbrush will work at next mode and start another 2 minutes cycle. The function turn of the button is: power on at normal mode, soft mode, massage mode and power off.

5. Apply light pressure and gently move the brush head slowly across the teeth in a small back and forth motion.

6. In order to maximize cleaning effectiveness, divide your mouth in 4 sections. Start with Section 1 and brush each section for 30 seconds. The brushing action will pause between each 30 second interval, indicating when to move on to the next section. The toothbrush will automatically stop after 2 minutes of use.(Fig. 11)

7. Use of the Interdental Brush Head

1). Specially designed shape is ideal for cleaning braces, crowns or spaces between poorly aligned teeth.

Fig.10


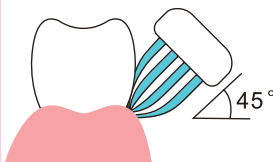
Fig.11

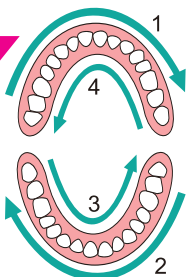2). Place the brush lightly and vertically against your teeth. Guide the brush head move slowly along the gap between the teeth.(Fig.12)

8. Pull brush head to detach from the tooth-brush.

9. Rinse both the brush head and handle under running water and wipe them dry.(Fig.13)

10. Put the handle back into the charger.

Fig.12


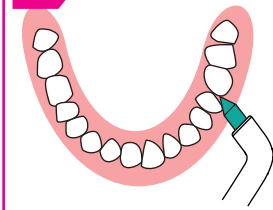
Fig.13

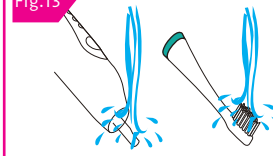
How to use Dental mirror

Use of the Dental Mirror

1. Allows vision in certain locations of the mouth where visibility is difficult or impossible.

2. By pressing the power button the LED mirror light will light up. Press again, the LED mirror light will go off.

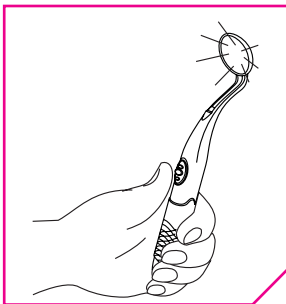Use only water or neutral cleaning agent to clean the device. Never use caustic or abrasive agents (e.g., vinegar or lime scale removers), since they may damage the device. Do not wash using hot water of over 40°C. Keep the device away from high temperature area and avoid direct sunlight.

Oral Irrigator

Main unit:

1. Remove the dirt by wiping using a cloth with liquid soap and wipe off using a dry cloth.
2. Do not submerge the main unit in water to clean.
3. Dry after washing with water.

Water tank cap:

Wipe off any grime and dirt with a cloth.

Nozzle:

1. Wash with water and wipe with a soft cloth.
2. Replace the nozzle every 6 months.

Water suction tube:

1. Wash with water and wipe with a soft cloth.
2. Do not bend, pull or twist the water suction tube.

Tank:

1. Wash with water and wipe with a soft cloth.
2. Dry the tank when not using for a long period.

Be sure to remove the plug from the electrical outlet before cleaning the product.

Sonic Toothbrush

Charger:

1. Pull the storage to remove.
2. Using a soft cloth, wipe off any dirt on the charger and the storage.
3. Clean the plug periodically to ensure it remains in good working order. The accumulation of dust or dirt on the plug can cause an electric shock, overheat or fire.

Toothbrush:

1. Remove the brush head.
2. Wipe all surface areas with a damp cloth and place the handle back into the charger. Do not submerge the handle into water.
3. Rinse the brush head and bristles after each use.
4. Do not clean in the dishwasher.
5. Any other servicing should be performed by an authorized service representative or that the product has no user serviceable parts.

ISO9001 & ISO13485 Certified Manufacturer

Designed By MTS in Iran, Assembled at China

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